

# injeTAK®

CYSTOSCOPY NEEDLE FOR PRECISE PLACEMENT



*Innovation for Pelvic Health*

# PATENTED DEPTH DIALING FOR ACCURACY, SAFETY AND SPEED

injeTAK®

## UNMATCHED PATIENT TREATMENT AND SAFETY

- For cystoscopy injections into the bladder wall for Overactive Bladder (OAB) and Neurogenic Detrusor Overactivity (NDO)
- Accurate, safe dispensing with decreased procedure time

## OPTIMIZED ADJUSTABILITY

- Retractable needle protects needle tip and cystoscope
- Does not require introducer sheath
- Needle tip can be adjusted to 0, 2, 3, 4 and 5 mm to accommodate bladder wall thickness
- Settings provide optimal tip length for speed, safety and consistency (2 mm per RX manufacture's PI)
- Ergonomic design allows for easy tip adjustments


## NEEDLE CHARACTERISTICS

- 4.8 Fr allows for use with smaller working channels and most pediatric scopes
- Black indicator tip provides visual proximity to tissue for accurate placement, safe penetration and gauged puncture depth
- Precise angle and 3-bevel tip engineered to minimize puncture force
- Unique cannula design provides easy injection and superior flexibility while maintaining structural integrity

## AVAILABLE SIZE OPTIONS

35 cm – DIS199

70 cm – DIS201

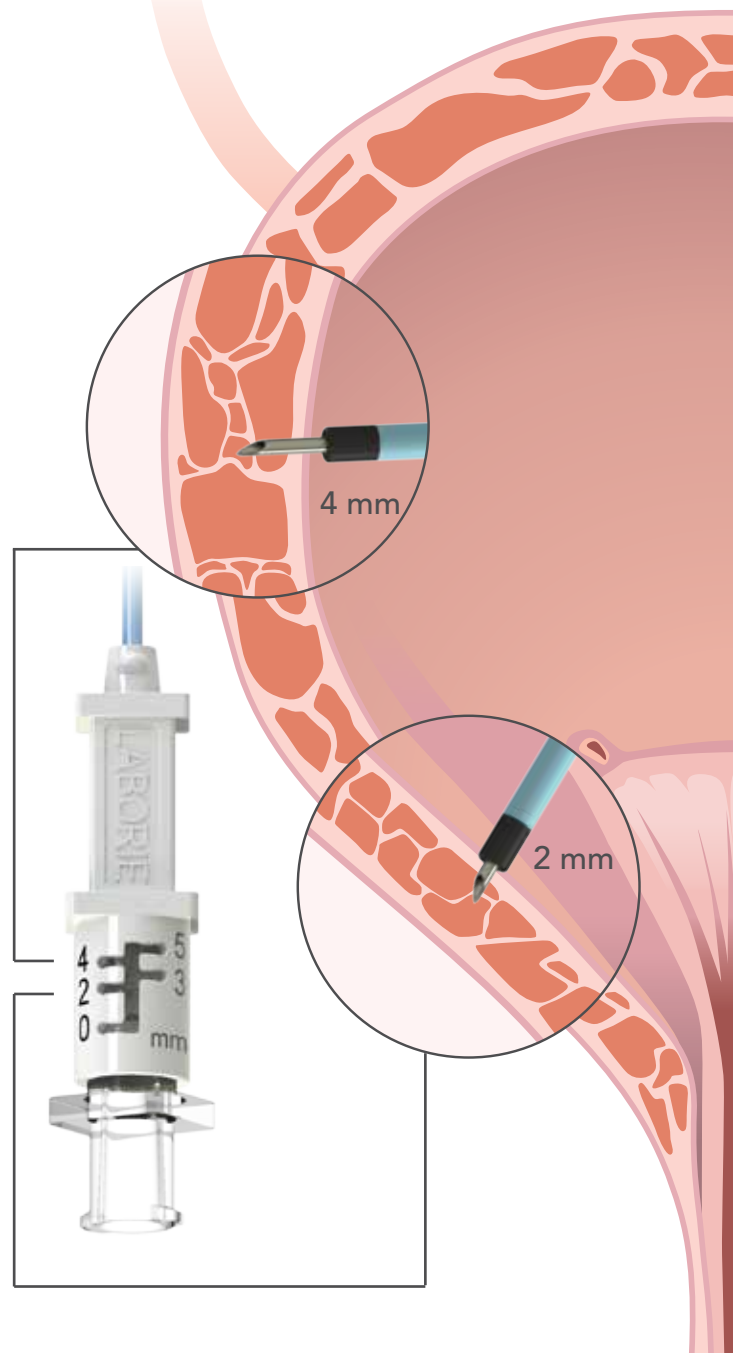

**USA:**  
Tel.: +1 802 857 1300  
800 522 6743  
Fax: +1 802 878 1122  
Email: [usmarketing@laborie.com](mailto:usmarketing@laborie.com)

**EUROPE:**  
Tel.: +44 1788 547888  
Fax: +44 1788 571652  
Email: [euromarketing@laborie.com](mailto:euromarketing@laborie.com)

**CANADA/INTERNATIONAL:**  
Tel.: +1 905 612 1170  
Fax: +1 905 612 9731  
Email: [intlmarketing@laborie.com](mailto:intlmarketing@laborie.com)

